Loo **Seal** **EVO**

DIY EDITION

User Manual

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Manuale d'uso

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Overview

User Manual (EN).....	3
Benutzerhandbuch (DE).....	17

Manuel d'utilisation (FR)
Manual del usuario (ES)
Manuale d'uso (IT)
Gebruikershandleiding (NL)
Podręcznik użytkownika (PL)
Uživatelská příručka (CZ)

Congratulations on your LooSeal EVO sealing toilet!

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Content (EN)

Special features of the DIY Edition	4
Dimensions of the toilet	4
Necessary distances.....	4
Weight distribution	5
Secure positioning - secure against slipping	5
1. Introduction	6
1.1 Packaging List.....	6
1.2 To this manual	6
2. Safety Instructions	7
2.1 Intended use.....	7
2.2 User qualification	7
2.3 Operation safety.....	7
2.4 Cleaning and maintenance.....	8
2.5 Modification and repairs	8
2.6 Residual risks	8
3. Design and function.....	9
3.1 Overview.....	9
3.2 Control panel & Display	10
4. Preparation	11
4.1 Check the condition of the toilet	11
4.2 Battery charging.....	11
4.3 Insert or replace liner	11
5. Operating steps after using the toilet	13
5.1 Use Super-Absorber	13
5.2 Note filling level.....	14
5.3 Start sealing process	14
6. Cleaning.....	15
7. Maintenance.....	15
8. Error Codes.....	15
9. Troubleshooting.....	16
10. Consumables and Accessories	16
11. Frequently Asked Questions (FAQ)	17
12. Disposal and Environment.....	17
12. Technical Specifications.....	18


Special features of the DIY Edition

Dimensions of the toilet

The body of the toilet has the following dimensions measured at the outermost elements:

Height: 24 cm (25 cm incl. silicone part)
Width: 36 cm (incl. pre-assembled leg axles)
Depth: 42 cm

Please note that additional buffer zones for opening the lid, for operating the control panel and for charging the toilet at the back must be taken into account for the functional use of the toilet. (see necessary distances)

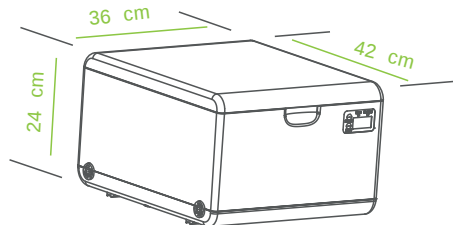Necessary distances

To ensure that you can continue to use the toilet with all its functions despite the individual construction, the necessary clearances must be maintained.

Distance from body to base:

There must be a minimum distance between the underside of the body and the floor to allow the flap to open fully and the contents of the bag to slide forward slightly. Without this distance and without the slightly slanted flap, the toilet will not function properly! The flap itself requires a distance of 12 cm until it is fully open. However, we recommend a total distance to the floor of at least 20 cm to prevent unwanted bag jams.

Distance at the rear:

Both the charging port and the main switch are located at the back of the toilet. The main switch should be accessible at all times so that the toilet can be switched off again after use. The battery life can only be used for the intended period of around 100 uses. Alternatively, the toilet can simply be charged using continuous power in standby mode.

The hollow plug supplied requires a little space for the charging port. When plugged in, the plug of the charging cable protrudes about 4 cm from the back wall of the toilet. Alternatively, an angled plug can be used instead of the supplied plug to reduce the space required.

Distance for opening the toilet lid:


The toilet lid is installed with a hinge that protrudes to the rear. If the toilet lid is opened by 90°, a buffer space of 3 cm must be planned behind the toilet. If the lid is to be opened fully, we recommend a buffer space of at least 4 cm. The lid requires approx. 37 cm of space at the top.

Clearance at the front of the toilet:

The control panel is located at the front of the toilet. If you want to conceal the front of the toilet, you should make sure that the toilet can still be operated via the display or the integrated buttons.


Weight distribution

With a self-built substructure, it is extremely important to load the toilet at the right points. There are 4 points at which the toilet may be loaded so that the toilet is not damaged and the maximum permissible load of 150 kg is not reduced. The weight load should be evenly distributed over these 4 points. You can recognize these 4 points by the pre-assembled silicone parts on the underside of the toilet (see illustration). All other points are not designed for loads and could break if installed incorrectly. The load could, for example, be applied via 2 rails/strips on the silicone components provided for this purpose.

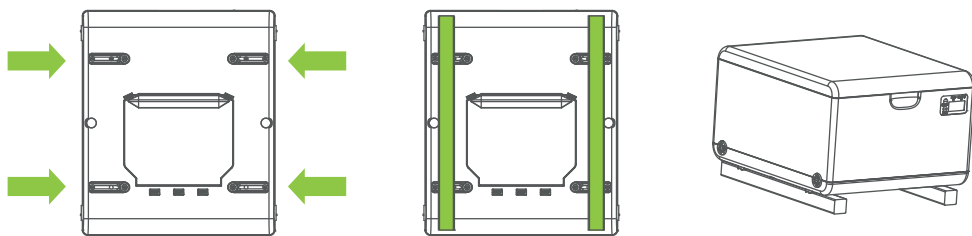Secure positioning - secure against slipping

The toilet should be secured against slipping during use and also in a moving vehicle. There are various ways to do this. We recommend using the screw openings in the silicone components and/or the openings on the sides to secure the toilet. Alternatively, the toilet can of course also be secured with an appropriate frame.

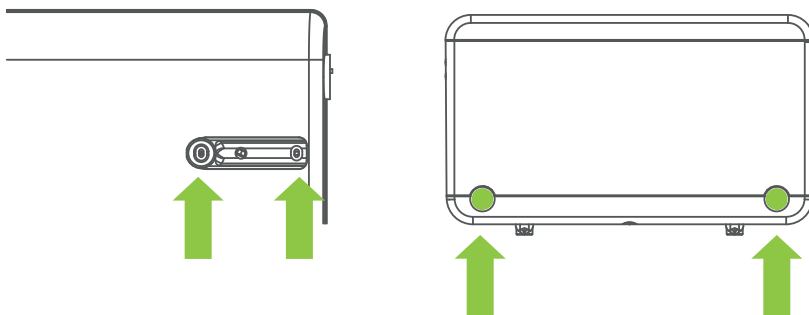
CAUTION:

The toilet body must not be drilled into or pierced with screws. There are often sensitive components behind the outer shell that can be damaged!

1. Introduction

1.1 Packaging List

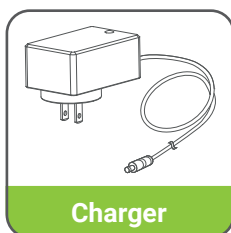Consumables:

1.2 To this manual

These operating instructions will familiarise you with the features and functions of the waterless sealing toilet LooSeal EVO. The instructions contain important information for the correct and safe use of the toilet.

- » Read these operating instructions carefully before using the toilet for the first time.
- » Keep these instructions in a safe place so that the information is available to you at all times.

These instructions are constantly being improved, so it is possible that the document enclosed with the LooSeal EVO may not correspond to the current version. We recommend that you check our website at <https://tomtur.com> to see if there is a newer version of this manual is available.

2. Safety Instructions

2.1 Intended use

The LooSeal EVO waterless sealing toilet is designed for the bagging of excrement and personal hygiene articles. The bagging of nappies or biowaste is also possible. Other uses are not permitted.

Possible misuse

Unintended use of the appliance may result in damage to property and personal injury. The appliance is not suitable for the following uses, for example:

- » Bagging of pointed or sharp-edged objects (e.g. shards, needles, razor blades)
- » Bagging of burning or glowing objects (e.g. ashes, cigarettes, matches)
- » Bagging of liquids and chemicals that could react with the bags
- » Bagging of food or living creatures
- » Use as a climbing aid

2.2 User qualification

Persons who are not familiar with the device, children and people with physical, sensory or mental impairments may only use the LooSeal EVO under supervision or under the instruction of a responsible person.

Repair work is reserved for authorised specialist personnel. Contact your dealer if necessary.

2.3 Operation safety

General:

- » Check regularly that all functions are being carried out correctly.
- » After each use and before leaving the LooSeal EVO unattended, make sure that the sealing process is completed and the toilet is ready for use.

The soft-close mechanism closes the toilet lid automatically. This mechanism is damaged if the lid is pressed down by hand.

- » Do not press the toilet lid down.
- » Never use the toilet without the foil liner inserted.

Built-in battery:

- » Charge the battery fully prior to initial use.
- » Only use the original charger approved by us to charge the battery.
- » DO NOT dismantle or modify the battery pack.
- » Switch off the power when toilet is lying idle and not in use.
- » It is recommended to charge the battery at least every 3 months.
- » Please stop using and switch off the toilet when battery is overheating, and call the dealer.

Correct use:

- » Toilet must be placed stably on a flat floor.
- » The Max load capacity is 150 kgs.
- » Keep the surroundings of the toilet free of objects, children and pets during toilet running.

- » Any parts of human body DO NOT enter the toilet, because moving parts may cause injuries.
- » DO NOT stand or sit on the toilet during toilet running.
- » DO NOT throw sharp-edged objects into the toilet.
- » DO NOT throw burning or glowing objects into the toilet.
- » DO NOT use the toilet in extreme conditions (e.g. extremely hot, cold, humid).
- » DO NOT use the toilet without filling or sealing one end of film.

2.4 Cleaning and maintenance

- » Switch off and disconnect the charger before cleaning and maintenance.
- » DO NOT clean the toilet with running water or immersing in water.
- » Use a soft moist cloth with neutral detergent to clean the toilet. DO NOT use bleach.
- » Cleaning parts: lid, seat, leg and surface of body. DO NOT clean any parts insides.
- » Maintenance shall only be carried out by qualified personnel authorized by the dealer or the manufacturer.

LooSeal foil liners ensure reliable functioning of your LooSeal EVO. The use of other foil liners can lead to malfunctions.

- » Use only original film liner and super absorber.

2.5 Modification and repairs

Repairs that are not carried out properly can result in accidents and damage to property.

- » DO NOT open the base body of the toilet.
- » If you suspect a defect, contact your dealer.
- » DO NOT install any additional components or modify the toilet.

2.6 Residual risks

Live parts

Live parts can be exposed by opening the base body. There is a danger of life-threatening electric shock.

- » DO NOT open the main body.

Water ingress

Water can penetrate the inside of the toilet through an unclosed lid or on through some toilet parts and damage the electronic control or the sealing mechanism. Malfunctions and component damage may result.

- » DO NOT expose the LooSeal EVO to jet water (e.g. when showering).
- » Only use the toilet while sitting down.
- » Never operate the LooSeal EVO in a wet room while humidity is very high.
- » If the sealing unit has come into contact with water: let the toilet dry for 24 hours without the foil liner inserted.

Substances hazardous to health




Ingestion of the absorber may cause nausea and stomach problems. In case of eye contact, there is a risk of severe eye irritation.

- » Keep the absorber out of reach of children and pets.
- » Observe the safety data sheet of the absorber, available on the website: <https://tomtur.com>
- » In case of eye contact, rinse the eye with running water for several minutes.

- » If the absorber is swallowed, rinse the mouth with water and drink plenty of water.
- » Consult a doctor if symptoms persist.

3. Design and function

3.1 Overview

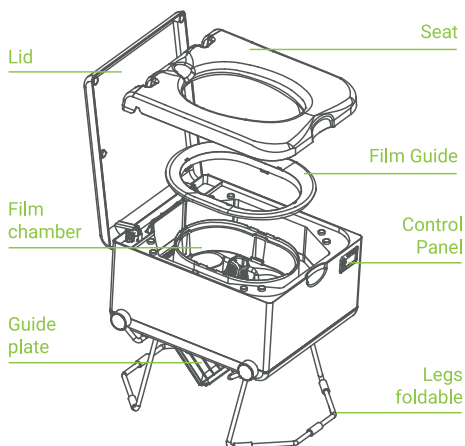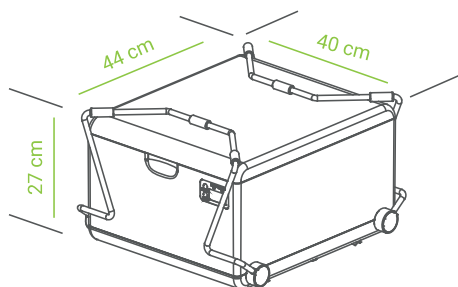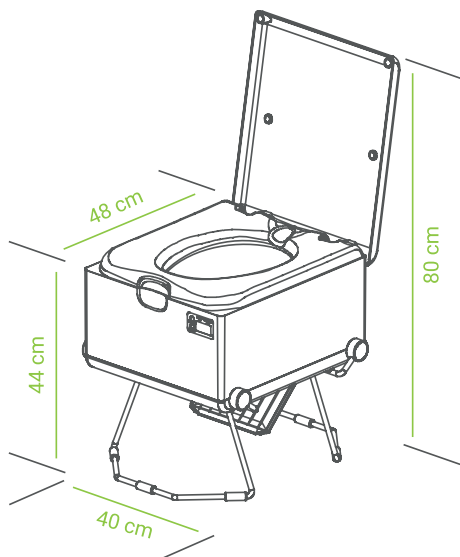The LooSeal EVO is the first portable sealing toilet that works on the principle of bag sealing. After each use, the toilet seals the toilet contents into a separate foil bag.

The LooSeal film liners are made from a continuous compostable foil, folded several times, which has previously been pulled onto the film guide and inserted into the film chamber. The foil pouch imitates a toilet bowl.


After triggering the sealing process via the „Start“ button on the control panel, the toilet pulls the liner from the foil guide automatically, seals it and prepares for the next use of the toilet.

After the sealing process, the odour-tight sealed foil bags fall out of the toilet at the bottom via the guide plate.

Key features:

- 1) Simple: one-key operation, even a child can operate it.
- 2) Hermetic: both ends of polymer bag are welded by heating to prevent from leaking.
- 3) Hygienic: the polymer bag wraps and seals all droppings and odors.
- 4) Independent: no need any connection to water supply and power supply.
- 5) Mobile: compact size and lightweight make it possible to be carried anywhere.
- 6) Sustainable: eco-friendly film is biodegradable, low energy consumption.

3.2 Control panel & Display

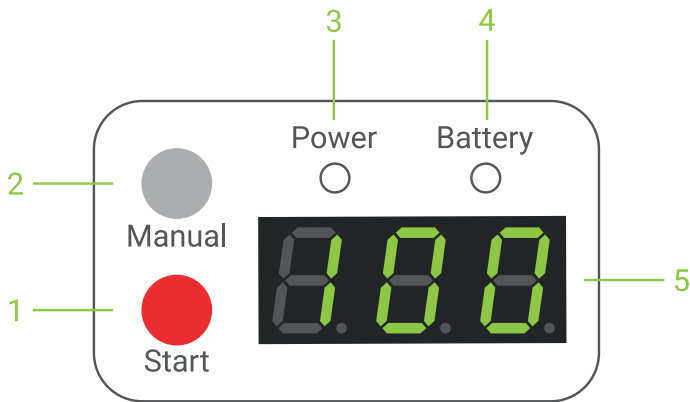No.	Description	Function
1	Button „Start“	Press and hold “Start” button for 2 seconds, and you will hear a beep sound, then the toilet starts sealing.
2	Button „Manual“	Press and hold the „Manual” button to pull the liner further down. Use case 1: To create larger bags. Use case 2: To change the liner.
3	LED „Power“	LED on means the power switch on the backside of the toilet is turned on. LED off means the power switch on the backside of the toilet is turned off.
4	LED „Battery“	LED Flashes green light: Battery is in charging mode. LED Steady green light: Battery is fully charged.
5	Display	Display in normal mode: The number on the display shows the remaining percentage of power in 25% steps. Operating mode / sealing mode: Countdown until the sealing process is completed. The system will enter sleep mode if there is no operation within 60 seconds, and can be woken up by pressing any buttons on the panel. (To save battery, it is recommended to turn the power switch on the backside of the toilet off if the toilet is not used.)

4. Preparation

4.1 Check the condition of the toilet

Before each use of the toilet, check if the toilet is in working condition.

1. The toilet is stable on both legs.
2. The guide plate is open. (This should happen automatically when the legs are positioned correctly)
3. The power LED is on or the power button on the back of the toilet has been pressed.
4. The battery indicator is above 25%. When the number is at 25% or lower, please charge it at soonest.
5. A liner is inserted correctly, tightly welded and not already at the end of its capacity. The end of capacity is recognized by the visible text imprint „STOP“ on the liner.
6. Make sure that there are always enough refill foil liners and LooSeal Super-Absorber near the LooSeal EVO toilet.

4.2 Battery charging

Use the original charger to charge the toilet (Input 100-240VAC, Output 14 VDC/2A). The power socket is located at the back of the toilet. The battery indicator on the control panel flashes while charging the battery, and it indicates steady green light when the battery is fully charged. The number on display means the remaining percentage of power. When the number is at 25%, please charge it at soonest. Please charge the battery before first use. The toilet can be used during charging if the battery capacity is high enough.

4.3 Insert or replace liner

The LooSeal liner is made of compostable film. This means that the liner begins to degenerate and leak a maximum of 6 months after opening the airtight package. We therefore recommend changing the liner after a maximum of 3 months, even if the liner has not been completely used up.


When you see“ STOP” on the film in the toilet, it means the film is going to run out, and you should refill a new film liner right now.

To replace an existing liner or to insert a new liner, please follow these steps:


1. Press and hold the “Manual” button on the control panel till the old liner in toilet fully runs out.


2. Take the seat and film guide out of the toilet. Place the film roll around liner base, and the arrow marks on the straps must point upwards.

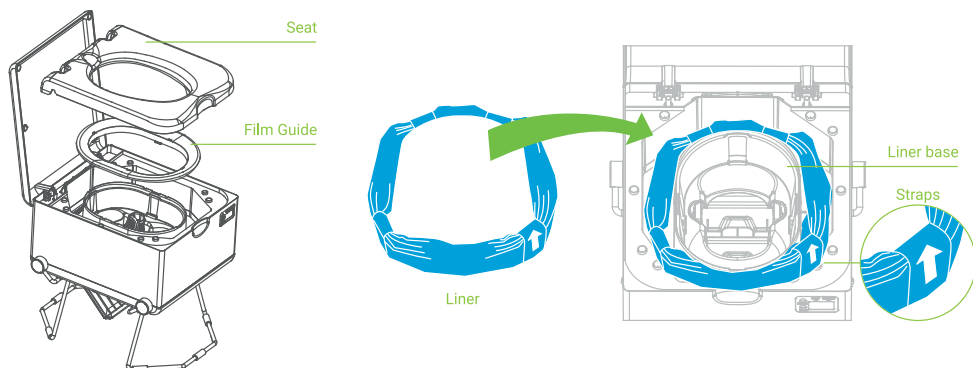
3. Remove the four straps, and place the film guide onto the liner base.

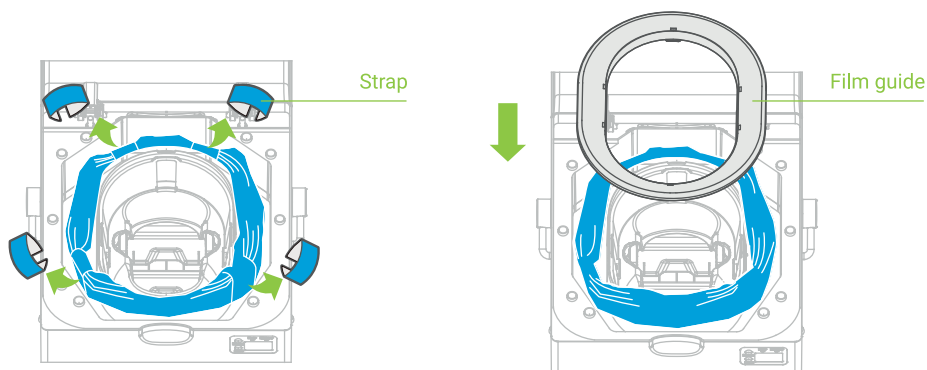
4. Pull the film upwards about 20cm from outer edge, then grasp it into a bunch and push it into the toilet. All the processes must be evenly done to avoid the liner misaligned feeding.

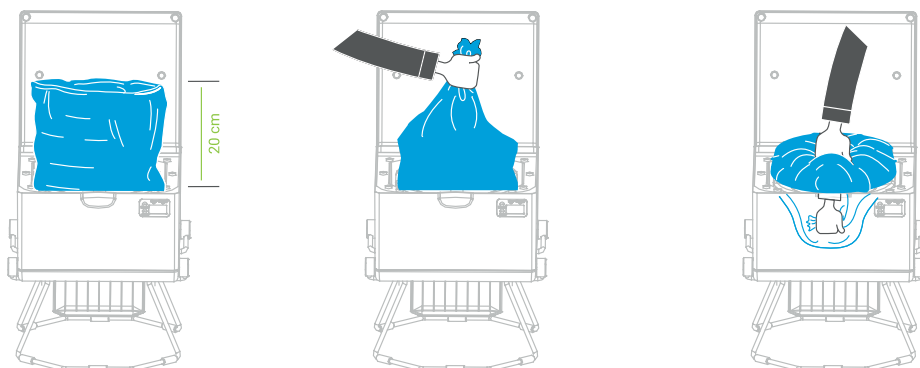5. Remove the jig from the holder at the rear end of the liner base. Align the jig with the arrow marks on the body, and push it with the film in between the two rollers at each side at the same time. Take out the jig and place it back in the indentation. Place the seat on the toilet body in correct position.

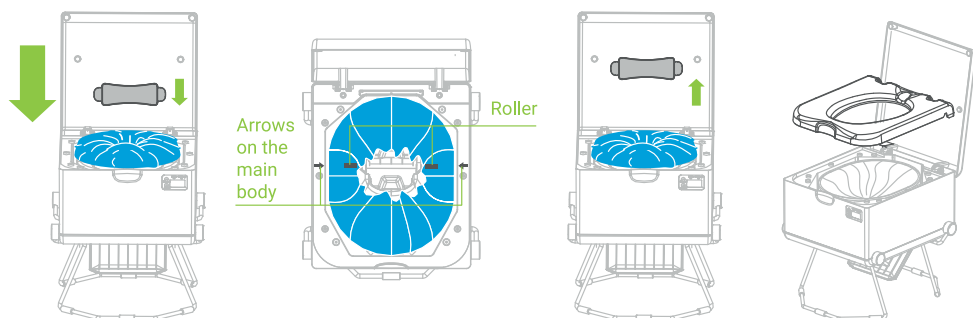6. Press and hold the “Start” button for 2 seconds, and you will hear a beep sound, then the toilet starts sealing. The numbers in display window are counted down from 100 to 00. You will hear a beep sound when it reaches 00, and the sealing is completed. The sealed bag or separated film will fall to floor. Take out the bag/film and dispose it like baby diaper according to local regulations and rules. The system will enter sleep mode if there is no operation in 60 seconds, and can be woken up by pressing any buttons on the panel.

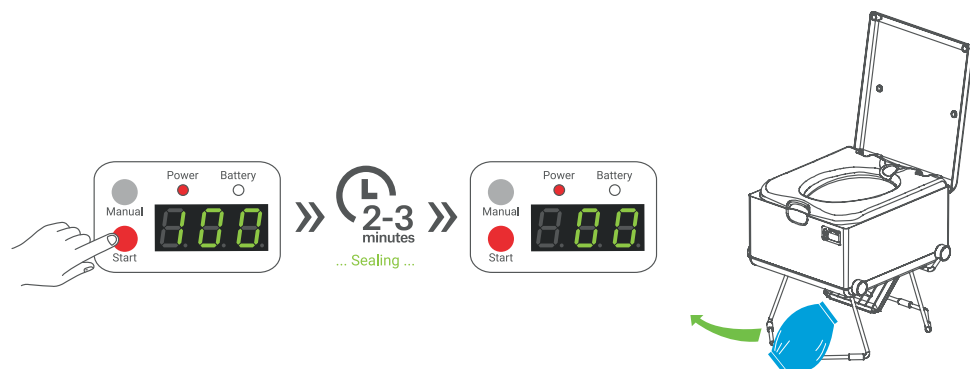
5. Operating steps after using the toilet

5.1 Use Super-Absorber



Add a bag of super absorber into the toilet before using the toilet. For binding liquids (e.g. urine) we recommend the use of LooSeal Super-Absorbers. The powder binds the liquids in the foil bag and thus serves as additional leakage protection, as well as odour prevention.

Super-Absorber

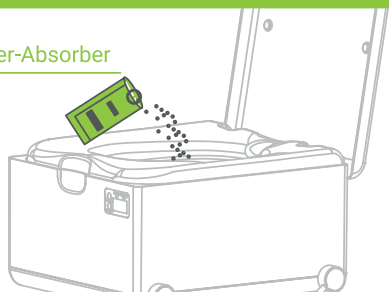
5.2 Note filling level

Toilet paper must not be above the specified maximum filling height, otherwise the weld seam will be interrupted. The maximum filling height is about half of the film pouch. If more space is needed in the bag or the bag needs to be enlarged, this can be done simply by pressing and holding the „Manual“ button.

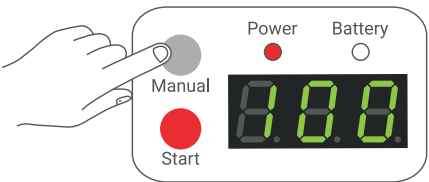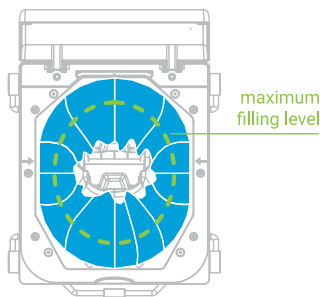5.3 Start sealing process

Press and hold the “Start” button for 2 seconds, and you will hear a beep sound, then the toilet starts sealing. The numbers in display window are counted down from 100 to 00. You will hear a beep sound when it reaches 00, and the sealing is completed. The sealed bag or separated film will fall to floor. Take out the bag/film and dispose it like baby diaper according to local regulations and rules. The system will enter sleep mode if there is no operation in 60 seconds, and can be woken up by pressing any buttons on the panel.

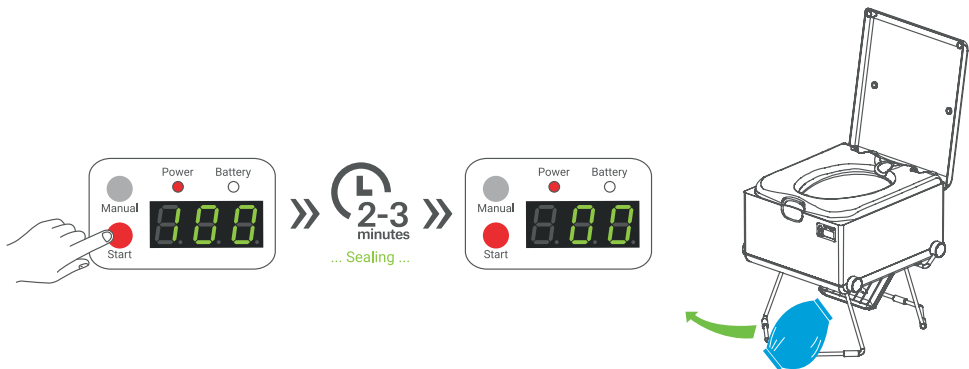6. Cleaning

As the LooSeal EVO works without water and the „foil bowl“ renews itself with every use, maintenance is very low. However, we recommend cleaning the seat, lid, leg and the surface of the toilet body regularly with a damp cloth and commercially available cleaning agents:

WARNING! Damage to the toilet may result from improper cleaning!

- 1) Switch off and disconnect the charger before cleaning and maintenance.
- 2) DO NOT clean the toilet with running water or immersing in water (e.g. with a hose).
- 3) DO NOT use bleach. Use a soft moist cloth with neutral detergent to clean the toilet.
- 4) DO NOT clean any parts insides.

7. Maintenance

Maintenance shall only be carried out by qualified personnel authorized by the distributor or the manufacturer.

WARNING! DO NOT open the base body.

Opening the base body can cause irreversible damage to the LooSeal EVO. This will invalidate the warranty.

8. Error Codes

Error	Fault	Solutions
E1	Thermosensor protection	Contact your dealer.
E2	Heating device protection	Contact your dealer.
E3	Timer protection There is foreign object in the clamping position.	Remove object and restart it by switching off and on. If the fault persists, please contact your dealer.
E4	Overcurrent protection when clamping	Remove object and restart it by switching off and on. If the fault persists, please contact your dealer.
E5	Overcurrent protection when declamping	Remove object and restart it by switching off and on. If the fault persists, please contact your dealer.


9. Troubleshooting

Fault	Possible Cause	Solutions
Control panel does not light	Power off	Press the power button on the back of the toilet.
	Low battery voltage	Charge the battery.
	Toilet in standby	Press any button.
	Control panel or toilet is defective	Contact your dealer.
Fail to charge the battery	Improper charger	Only use the original charger.
Power off while sealing	Low battery voltage	Charge the battery.
Unstable	Legs are not fully open	Open legs till you hear a click.
Separation not successful	There is a foreign object in the clamping position	Remove the object and switch the toilet off and on. Then start the sealing process again.
	Heating or cutting device is defective.	Contact your dealer.
Broken film	-	Press and hold „Manual“ button to feed the film till no more holes on in the film. Press the „Start“ button to seal the film as described in the chapter „Preparation“.
The film wraps itself around the rollers	The bags accumulate on the floor and are pulled upwards and in the rollers.	Remove the film from the rollers. Remove the bags on floor. Press the „Start“ button to seal the film as described in the chapter „Preparation“.

10. Consumables and Accessories

Only use original products on the toilet. Original consumables and accessories can also be purchased quickly and easily from the ToMTuR online shop:

<https://tomtur.com>

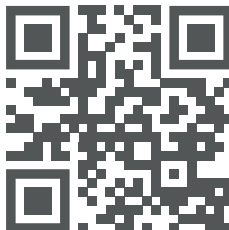11. Frequently Asked Questions (FAQ)

You can find all FAQ & information on <https://tomtur.com>

Can the bags be put on the compost (organic waste)?

Although the liner film is compostable, the bags should nevertheless be disposed of in the residual waste. Please pay attention to the regional laws that also apply to e.g. diapers.

Are the bags odour-proof?

At room temperature, the bags are odor-proof for between a few hours and up to 5 days, depending on the external conditions and contents. At higher temperatures, it is advisable to dispose of the bags within 2 days. We generally recommend disposing of the bags as soon as possible. The odor tightness has no influence on the liquid tightness of several weeks.

May I use the toilet again immediately after the sealing process has been started?

No, the previous „sealing process“ must be completed. Otherwise, malfunctions or defects may occur.

Can I also seal other things in the bags?

Yes, in principle organic waste, nappies or feminine hygiene products can be sealed in. However, no sharp-edged objects (broken glass, etc.) or hot ashes may be placed in the bags. Observe the maximum filling height!

Does the toilet work with a voltage of 12 volts?

The toilet works with a lithium-ion battery with a 12V base, but is charged by a 14V, 2A power supply (28W). In principle, it is possible to permanently connect the toilet to an existing 12V supply. Please note the maximum power supply of 28W!

How stable is the toilet?

The load capacity of the toilet is 150 kg.

12. Disposal and Environment

The LooSeal EVO complies with the requirements of the EU Directive „Restriction of Hazardous Substances“ (2011/65/EU). It is largely free of environmentally hazardous substances such as lead, cadmium, mercury or chromium VI.

In accordance with the EU Directive, the appliance is classified as waste electrical and electronic equipment electronic waste and must not be disposed of as household waste.

- Dispose of the product in accordance with local regulations.
- Dispose of used batteries at the designated collection points.

12. Technical Specifications

Attribute	Value
Model	VT-LS-001-DIYG
Colour	Grey
Dimensions	42 x 36 x 24 cm (DxWxH)
Recom. room width	40 cm
Recom. seat height	44 cm
Net Weight	7 kg
Max. Load Capacity	150 kg
Material	ABS, PP, POM
Battery	Lithium-ion 57 Wh
Battery Charging Time	approx. 2 h
Battery Charger	Input 100-240 VAC, Output 14 VDC/2 A
Battery Range	approx. 100 uses (1-2 weeks with 2 persons in permanent use)
Power	60 W (max)
Power consumption (per toilet use)	0,4 Wh
Power consumption (Standby)	0,3 Wh
Operation condition	0°C - 40°C, 20% - 80% RH
Storage conditions	-20°C - 60°C, 5% - 95% RH
Waterproof	No
Handles	Yes
Toilet Lid	with soft close
Cycle for Sealing Process	approx. 2-3 minutes
Recom. Substrate	LooSeal Super-Absorber
Mounting	Owners choice
Disposal of Bags	Residual waste

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ToMTuR GmbH